

QUICK CARD

Clinical Posture Screening

How it works – one page

Clinical Posture Screening is the TrainersEye posture workspace, run by the MMSx Posture engine.


TrainersEye app → Dashboard → Clinical Posture Screening · Clinical plan

GETTING STARTED

INSTANT SCAN

DETAILED ASSESSMENT


REPORT AND PRINT

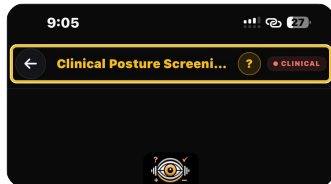

1

Open the tool


From the Dashboard, tap Clinical Posture Screening.

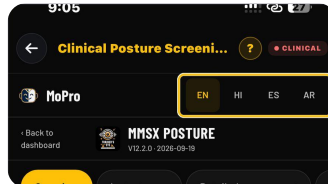

2

Choose your language


Pick EN, HI, ES or AR beside the MoPro header.

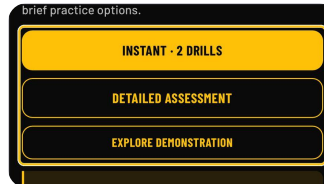

3

Three ways in


The Overview explains the method – see the visible differences first, add another view when useful.

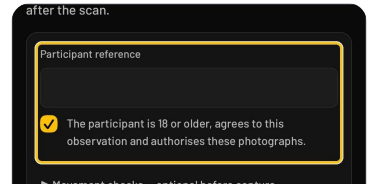

4

Consent first


Instant scan starts with a photograph, but consent comes first.

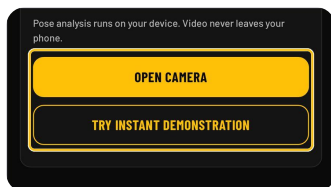

5

Capture the views


Tap Open camera and photograph the participant from the front, the right side, the left side and the back.

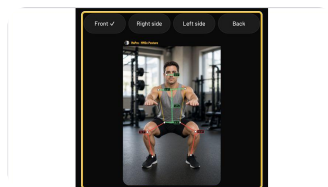

6

Read a view


Each captured view shows the pose lines drawn on the photograph – ear line, shoulder line, trunk shift.

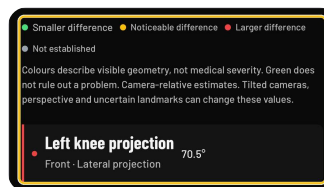

7

Read the findings


Under the photograph, every finding is a tile: its name, the view and projection it comes from.

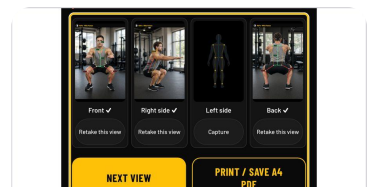

8

Your observation summary


The summary shows the four views with a tick on each one captured, Retake this view under each and Capture under any still missing.

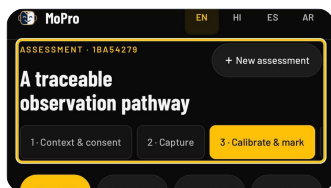

9

Start the traceable pathway


Detailed assessment records an assessment with its own ID and walks a numbered pathway: 1 · Context & consent, 2.

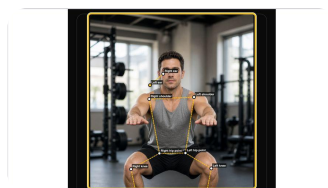

10

Calibrate and mark


In Calibrate & mark, each view carries named landmarks – right and left ear, shoulder, hip point, knee and ankle.

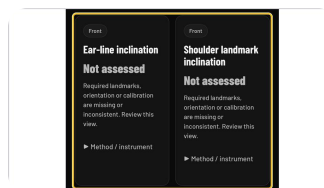

11

The measurement cards


Every measurement is a card by view – Ear-line inclination, Shoulder landmark inclination.

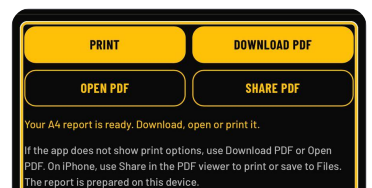

12

Print or save the report

Print / Save A4 PDF opens the A4 report preview with four buttons: Print, Download PDF, Open PDF and Share PDF.

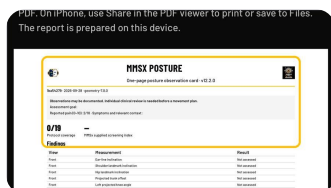

13

The one-page observation card

The detailed assessment prints its own card: protocol coverage out of nineteen and the MMSx screening index.



You're ready

That's it. Open Clinical Posture Screening and capture the first view.

ON THE DOCUMENT

REPORT MASTHEAD · RECORD LINE · STATUS BANNER · FOUR VIEWS · VISIBLE POSTURE FINDINGS · LEGEND · GEOMETRY NOTE · MOVEMENT PROGRAMME · METHOD LINE · NOT ESTABLISHED